MOBILE LIFTER FOR ATEX ENVIRONMENTS

Max 300 kg



TABLE OF CONTENT



Svenema develop, design and manufactures custom-made ergonomic lifters.

Our focus is to manufacture mobile lifters with an ergonomic design for material handling in the pharmaceutical environment.

We have provided the pharmaceutical industry with mobile lifters since 1995.

ABOUT ATEX LIFTERS Can be used in zone 1, 2, 21 and 22.

EXAMPLE OF MODULES 6 Inspirational pictures of what goods you can handle with the mobile lifters for ATEX environments.



TECHNICAL DATA Max 300 kg. Max lifting height 2000 mm.


WHY CHOOSE SVENEMA? Always in stainless steel 304 or 316. Ra finish 1,4 or better.

GALLERY Some of the custom-made mobile lifters we've done through the years.

Svenema AB • info@svenema.se • www.svenema.se





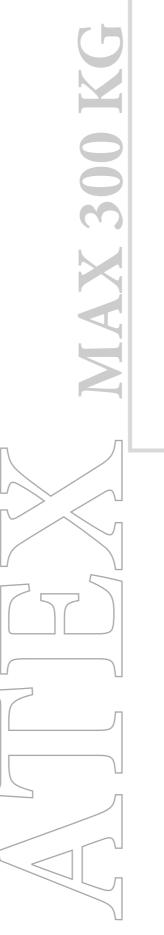



ABOUT ATEX LIFTERS

Our ATEX rated lifters can be used in zone 1, 2, 21 and 22.

Instead of a lifter driven by electricity the ATEX lifters are driven manually or with a pneumatic solution.


All our ATEX lifters can be customized with the same functions as our regular mobile lifters, such as tilting function, clamp function, integrated scale, etc.





Handling of filter


- Manual lift by winch and crank.
- Manual clamp to hold the filter.

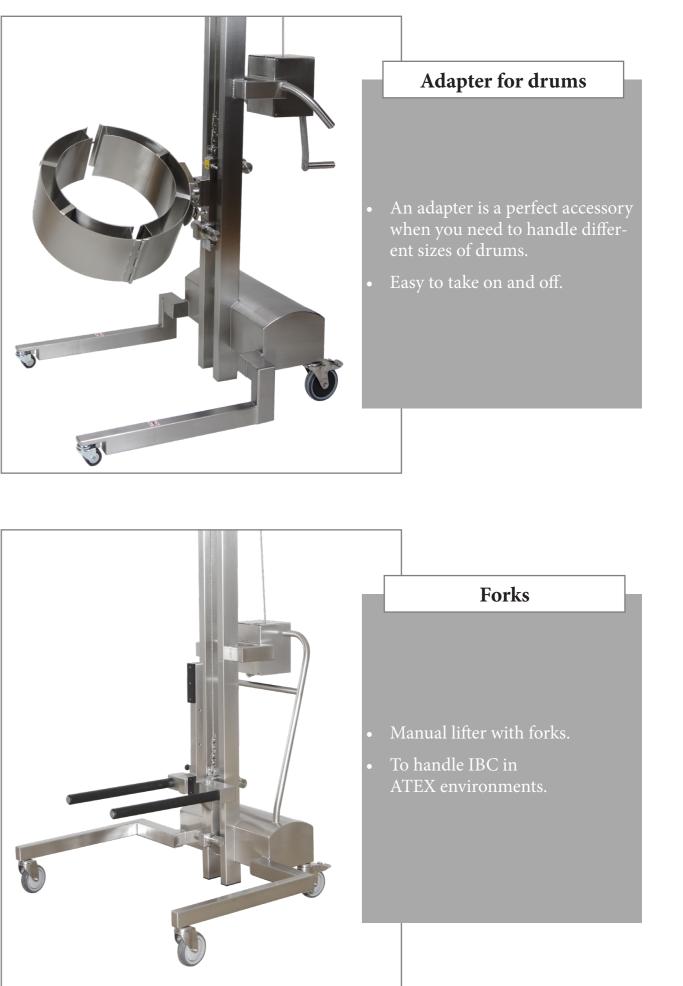



Handling of drum

- Manual clamp to hold the drum.
- winch and crank.

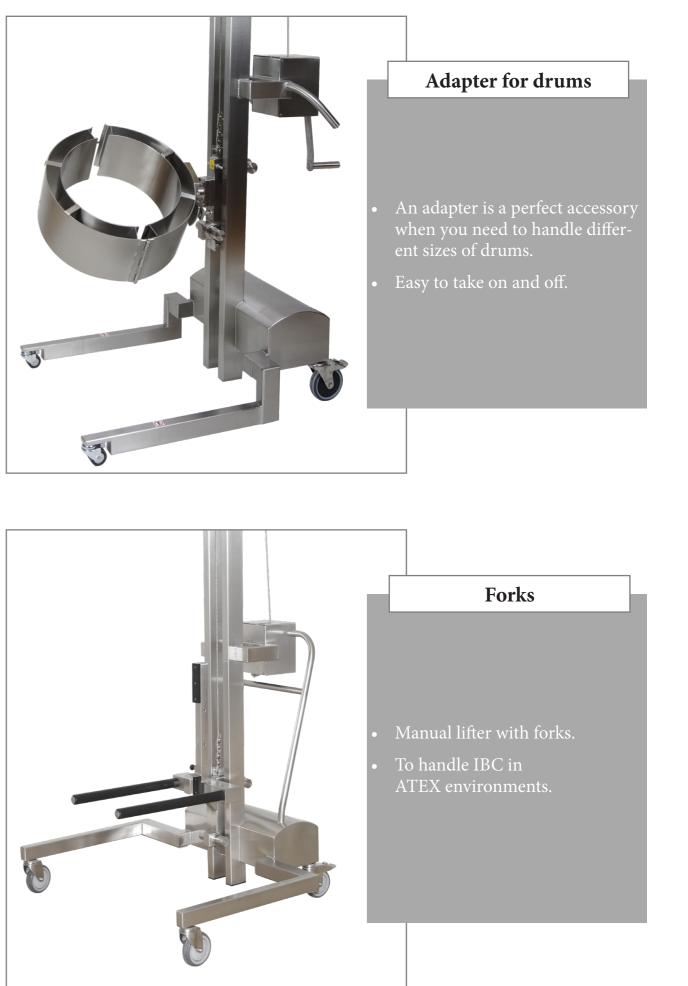




ATEX column lifter

Pneumatically driven.

- With a pneumatic column lifter you have the ability to lift good very heigh.
- Ideal for lifting goods to a mezzanine level.

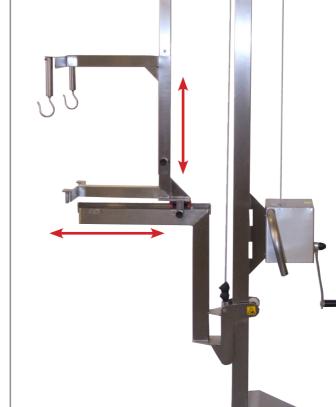




Special designed arm

- To handle drums from pallet to a scale in an isolator.
- Telescopic arm to reach into the isolator.
- The hoop holds the drum in place, making sure it can't wobble.



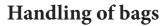

- To handle bags of different sizes.
- The hooks are flexible so that the bag will always be in stretched position even when it gets more and more empty.
- The telescopic movement helps you to reach higher and longer when needed.



- A small and compact lifter that handles drums.
- Lift the drum straight up and move it from one place to another.
- Manually driven by a winch.

TECHNICAL DATA

- Max 300 kg.
- Max lifting height 2000 mm.
- All in stainless steel 304 or better. •
- Stainless steel 316 with certificate • is optional.
- Surface finish Ra 1.4 or better. Ra 0.8 is optional.
- Smooth and clean welding, grinded and polished as an option.
- Smooth surfaces no entrapment areas.

- All components are industrial grade components.
- ATEX zone 1, 2, 21 & 22.
- Manual or pneumatic solutions.
- Must be connected to ground or used on conductive floor.





WHY CHOOSE SVENEMA?

- Always in stainless steel 304 or 316.
- Ra finish 0,8 1,4.
- Hygienic design easy to clean, minimizes the risk of cross contamination.
- cGMP.

- Specialized in cleanroom environments in the pharmaceutical industry.
- Only industrial components.
- Very low maintenance.
- Many years of experience of the pharmaceutical industries demanding environments.







Pneumatically driven column lifter. Can be used in ATEX zone 1,2,21 and 22.

Column lifter with electrical lift and tilt. Can handle containers up to 300 kg. The column lifter can be bolted in the floor or on the wall.





Our heavy duty lifter MAXI 300 is a powerful and reliable lifter. With its robust design it can lift goods up to 300 kg.

Mobile lifter midi to handle sterile components into a syringe filling line.

GALLERY

Some of the custom-made lifters we've designed through the years.

Mobile lifter for ATEX zone 1,2,21 and 22. We offer manually driven lifters or pneumatic solutions. No source of ignition.



Manual clamp and tilt for handling of drums. With the help of an adapter you can lift different sizes of drums.

SVENEMA AB PHARMACEUTICAL EQUIPMENT SINCE 1995.

